



Quick Start Guide

Pressio[®] 2 System

ICP/ICT Monitoring





Summary

Monitor set up and powering on	3
Catheter zeroing	4
Patient Monitor calibration	5
Settings	6

 Please refer to the Instructions For Use - NT500

○ Monitor set up and powering on

- 1 Affix Pressio® 2 Monitor to a pole



- 2 Connect Pressio® 2 Monitor to power


- 3 Connect Catheter Extension Cable to Catheter outlet (white)



- 4 Power on Pressio® 2 Monitor by pressing the [ON] button



- 5 Pressio® 2 Monitor performs a 4 seconds self-test

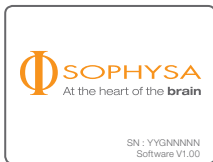

Catheter zeroing




Zero the catheter **BEFORE** Implantation

Pressio® 2 Monitor displays :

Connect catheter

- 1 Prepare a shallow cup with sterile saline within a sterile field

- 2 Connect Pressio® Catheter to the Catheter Extension Cable (up to the blue line)

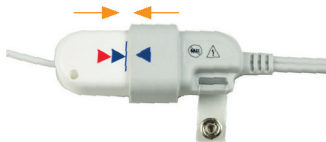

- 3 Immerse tip of Pressio® Catheter into sterile saline, without touching cup



- 4 Press [Zero] button  on Pressio® 2 Monitor touchscreen

- 5 Wait 4 seconds while Pressio® Catheter completes zeroing

Catheter is zeroing

70%

Zero catheter OK
Implant catheter



Patient monitor calibration

- 1 Connect Pressio® 2 Monitor to **patient** monitor with the Pressure Cable and Temperature Cable



2 Sending signal
0 mmHg
to patient monitor

- 2.1 Press [Zero] button on the **patient** monitor

- 2.2 Confirm **patient** monitor displays a steady 0 mmHg

- 2.3 Press [Ok] on Pressio® 2 Monitor touchscreen



3 Sending signal
30 mmHg
to patient monitor

- 3.1 Confirm **patient** monitor displays a steady 30 mmHg

- 3.2 Press [Ok] on Pressio® 2 Monitor touchscreen



4 Sending signal
60 mmHg
to patient monitor

- 4.1 Confirm **patient** monitor displays a steady 60 mmHg

- 4.2 Press [Ok] on Pressio® 2 Monitor touchscreen




Settings - Language




- 1 Press [Menu] button



- 2 Access the language menu by pressing [Language]

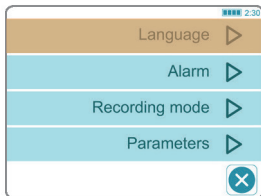

- 3 Select the language by pressing on the desired language

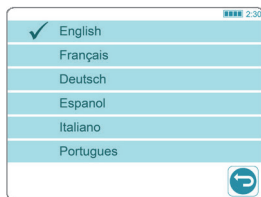

- 4 Press [BACK] to return to menu



- 5 Press [EXIT] to return to monitoring




Settings - Alarms

Activate / Disable | Thresholds

1 Press [Menu] button



2 Access alarm menu by pressing [Alarm]

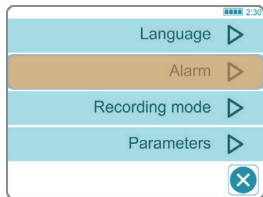


3 Press [ON] to activate alarm



3.1 Press [OFF] to disable alarm (PIN 1234)



3.2 Adjust the thresholds by pressing [+] or [-]

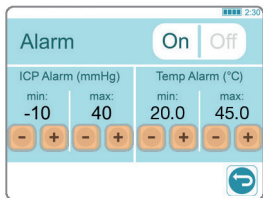

4 Press [BACK] to return to menu



5 Press [EXIT] to return to monitoring





Sophysa SA :

5, rue Guy Moquet
91400 Orsay - France
Tel.: +33 (0)1 69 35 35 00
Fax: +33 (0)1 69 35 36 90
contact@sophysa.com

Sophysa SA :

Chausée Bara 68
1429 Braine-l'Alleud - Belgique
Tel.: +32 (0)23 87 19 48
Fax: +33 (0)23 87 46 83
infobenelux@sophysa.com

Sophysa USA :

503 E Summit Street, Suite 5
Crown Point, IN 46307 - USA
Tel: +1 219 663 7711
Fax: +1 219 663 7741
contactusa@sophysa.com